


Sperm collection

SPERMACUP

Contenitore per raccolta
liquido seminale umano
Human semen collector

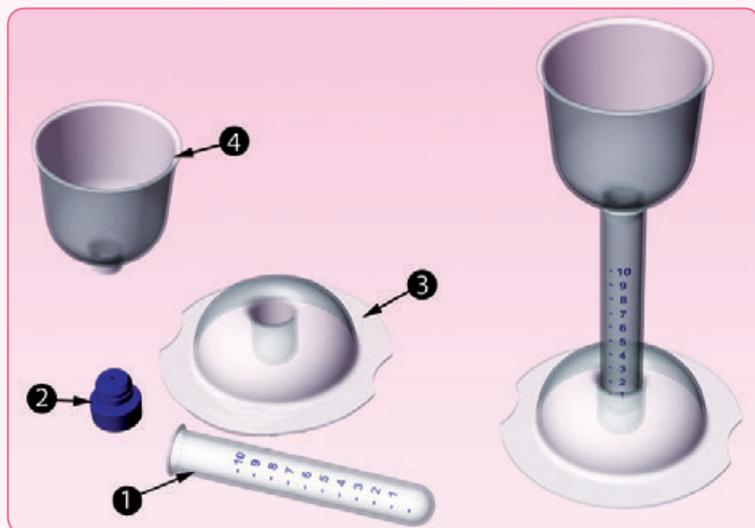

► Permette una raccolta di sperma umano senza manipolazione. Il dispositivo è composto da:

- una base **(3)**;
- una provetta graduata da 10 ml **(1)** con relativo tappo **(2)**;
- una coppa di raccolta **(4)**.

Montaggio semplice e veloce: inserire la provetta stabilmente nella base e la coppa di raccolta a sua volta nella provetta. Il dispositivo permette di eliminare la procedura di travaso dello sperma che in questo caso viene raccolto direttamente nella provetta. Sterilizzato a Raggi Beta per escludere ogni rischio di tossicità causato da eventuali residui di Ossido di Etilene.

✚ It allows a safe and direct human sperm collection without handling. The device is composed by:

- a basement **(3)**;
- a 10 ml graduated test tube **(1)** with plug **(2)**;
- a collection cup **(4)**.

Easy and fast assembly: first insert the tube firmly into the basement and then the collection cup into the tube. The device allows direct collection into the test tube and avoids sperm transfer.

Beta irradiated to eliminate any risk of toxicity from residual contamination common with EtO.

Cod.

770213

Scatola da 50 pezzi.

Fornito in confezione singola, monouso, sterilizzato a raggi Beta e marcato CE.

50 pcs. per box.

Single packaged, single use,
Beta sterilized and CE marked.